

Bioresonance Mk 6 – XXX Pro



Features

30 Watt Bioresonance Output.

1 to 10,000 Hz frequency adjustment .(100 to 2000 Hz on Mk 5)

594 Pre-programmed treatments

100 User- programmable treatments

210mm wide(without handle),105mm high,and 280mm deep.

Power : 1A @ 230V a.c. , 2A @ 110V a.c.

External fuse : 2A for 230V a.c. , 3.15A for 110V a.c.

LCD display light.

Removable handle for carrying or tilting on bench.

Case vents and micro fan for cooling.

Beeper for pushswitch presses and acupuncture sensor.

Adjustable "Deviation" control for programmable frequency sweep across centre frequencies.

Accessories(all can be used with Mk6)

Laser unit - 3mW ,670nm red laser for skin and acupuncture point treatments.(Green also now available)

Bioresonance modulated coloured light units for general bathing of tissue with Bioresonance light.(red,green,yellow,orange,blue, Infra-Red)

Energiser unit for Bioresonance energised water and homeopathic solutions.

Aqua unit for Bioresonance energised bathing liquids (immersion)

Micro Applicator for small or difficult to reach areas.

Acupuncture sensor / electro treatment unit (available soon).

Cymaresonance Ltd

14 The Moorlands , Malvern, Worcestershire , UK

Telephone : 361 288 4502 (USA)


+44 121 288 1022 (UK)

+44 7778 913102 (UK mobile)

Website : www.bioresonance.co.uk

email : sales@cymaresonance.co.uk



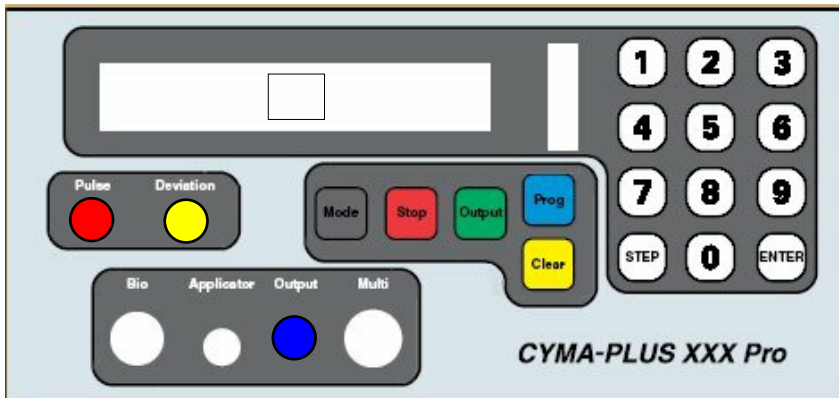

Bioresonance Mk 6 - XXX Pro Instructions

Power on

Press push button power switch on rear panel.

Green and red LED arrays should pulse up and down and lowest green LED segment should stay lit.

Two short beeps should sound to show internal power circuits are OK. If beeps carry on, there is a fault and the unit should be turned off and returned for investigation.

Message LCD display should show "CYMATICS Mk 6" briefly before showing "Program?"

Programs

Enter program number followed by "Enter" (#)

The LCD display will show the program number and channel (user programs 1-99) or condition (pre-set programs 100-594).

On user programs, press the yellow "Clear" button to clear the channel before entering frequency (up to 9999Hz and then "Enter" (#).

Press "Step" () to change to next channel of that program.*

On pre-programmed programs, press "Step" to go up one program.

To change program number, press the blue "Program button.

Outputting

To output the selected program, press the green "Output" button.

The LCD display will show the Program number and "Output>".

To stop output, press the red "Stop" button.

The LCD will then show the Program number and channel (user prog) or condition (pre-prog.)

Any button press will cause a short beep. If not, there is a fault.

Adjust output using the blue "Output" control.

The red LED array will show output level.

Adjust pulse duration if required using red "Pulse" control.

The yellow "Deviation" control gives an programmable regular frequency sweep on either side of the set frequencies.

The large Applicator, small Applicator, Energiser, Chroma-light units and Aqua units all use the 7mm jack socket. (second from left).

The foot plates use the two pin socket on the left.

The lasers use the 7 pin DIN connector (3rd from left).

Voltage setting for 120V/240V is by means of an internal switch which is accessible if the 4 screws holding the top lid are removed and the lid taken off.

USE OF THE UNIT AT THE WRONG VOLTAGE WILL RESULT IN DAMAGE

Use Earthed supplies for safety.

Do not obstruct the air intake grill on the rear of case or louvered grills on top of case, or overheating of the Amplifier may occur at high output levels.

Turn output off before switching off.

If unusual characters appear on the LCD display and/or if the machine "locks up", or if you lose some of your user programs after switching the machine off and on, it may be that the internal battery is flat (Approx. 5yrs) or that the memory has become corrupted due to power line surges.

The user memory can be cleared by pressing the yellow "Clear" button immediately after

Power up. The display should come up with "Clear RAM?". Press "Enter" (#) to confirm this.

The display will change to "Clearing RAM" for about 15secs before going to "Program?"

This reset will erase any user programs.

Manual Scan

It is possible to alter individual channel frequencies actively during output to experimentally fine tune commutations.

This feature is available during output of the 1-99 user programs but not the pre-programmed 100 upwards programs.

Press green Output button to start output and then press grey Mode button.

The LCD will show "Manual Scan?"

Press "#" and the LCD will show one of the 5 channel frequencies for that user program.

To increase the frequency by 1Hz, press the blue Prog. Button.

To decrease the frequency by 1Hz, press the yellow Clear button.

You will hear the harmonic frequency change in the output applicator.

To go to another channel, press the "*" (Step) button.

Press the grey Mode button to go back into normal Output mode.

***These alterations to your user program will be saved automatically to that program number.**

Active Program Step

Another feature you may use on the 1-100 user programs is Program Step where you may go up or down a program actively.

This could be useful if you need to go to another user program without interrupting a treatment.

Press green Output button to start output and then press grey Mode button twice.

The LCD will show "Go to next Program?"

Press "#" and the LCD will show one of the 5 channel frequencies for that user program.

To increase the Program number by 1, press the blue Prog. Button.

To decrease the Program number by 1, press the yellow Clear button.

You will hear the harmonic frequency change in the output applicator.

Press the grey Mode button to go back into normal Output mode.

Both Manual Scan and Active Program Step can be used during the same treatment output provided you press the grey Mode button to go back to normal output.

Pre-Programmed Commutations

100	ABSCESS	154	BACH FLOWER CERATO
101	ACARUS SCABIEL	155	BACH FL. CHERRY PLUM
102	ACIDOSIS	156	BACH FLOWER CHESTNUT
103	ACIDOSIS ACID/BASE	157	BACH FLOWER CHICORY
104	ACNE	158	BACH FLOWER CLEMATIS
105	ACONITE	159	BACH FL. CRAB APPLE
106	ACUPUNCT.BLOCK(REM)	160	BACH FLOWER ELM
107	ACUTE HEARING	161	BACH FLOWER GENTIAN
108	ACUTE SENSE OF SMELL	162	BACH FLOWER GORSE
109	ADDISONS DISEASE	163	BACH FLOWER HEATHER
110	ADENOIDS	164	BACH FLOWER HOLLY
111	ADHESIONS	165	BACH FL. HONEYSUCKLE
112	ADHESIONS(ALT)	166	BACH FLOWER HORNBEAM
113	AIR	167	BACH FL. IMPATIENS
114	ALLERGY	168	BACH FLOWER LARCH
115	ALOPOECIA	169	BACH FLOWER MIMULUS
116	ANAEMIA	170	BACH FLOWER MUSTARD
117	ANAEMIA(PERNICIOUS)	171	BACH FLOWER OAK
118	ANGINA(A)	172	BACH FLOWER OLIVE
119	ANGINA(B)	173	BACH FLOWER PINE
120	ANGINA PECTORIS	174	BACH FL.RED CHESTNUT
121	ANKYLOSIS(A)	175	BACH FL. ROCK ROSE
122	ANKYLOSIS(B)	176	BACH FL. ROCK WATER
123	AP MAG	177	BACH FL. SCLERANTHUS
124	APONEUROSES(MUSCLES)	178	BACH FL.STAR OF BETH
125	APPENDICITIS	179	BACH FL.SWEET CHESTN
126	ARACHNOID MATER VILL	180	BACH FLOWER VERVAIN
127	ARTERIES(ATHEROSCLE)	181	BACH FLOWER VINE
128	ARTERIO SCLEROSIS(A)	182	BACH FLOWER WALNUT
129	ARTERIO SCLEROSIS(B)	183	BACH FL.WATER VIOLET
130	ARTHRITIS-GENERAL	184	BACH FL.WHITE CHESTN
131	ARTHRITIS-DEGEN-OSTE	185	BACH FLOWER WILD OAT
132	ARTHRITIS-OSTEO	186	BACH FL. WILD ROSE
133	ARTHRITIS-DEGEN	187	BACH FLOWER WILLOW
134	ARTHRITIS-RHEUMATOID	188	BACH FL.RESCUE REMED
135	ARTHRITIS-ACUTE	189	BALANCE (VERTIGO)
136	ARTHRITIS-ELBOW	190	BELLS PALSY
137	ARTHRITIS-SACRUM	191	BENIGN TUMOR
138	ARTHRITIS-KNEES	192	BETA ENDORPHINES
139	ASTHMA	193	BILE(GALL BLAD.COND)
140	ASTHMA-BRONCHIAL(A)	194	BLADDER
141	ASTHMA-BRONCHIAL(B)	195	BLADDER IRRITATION
142	ASTHMA-CARDIAC(A)	196	BLISTER
143	ASTHMA-CARDIAC(B)	197	BLOCKAGE
144	ASTHMA-HAY(A)	198	BLOCKED GLAND
145	ASTHMA-HAY(B)	199	BLOOD PRES.-HIGH(A)
146	ASTHMA-BREATH SPASM	200	BLOOD PRES.-HIGH(B)
147	ATAXIA	201	BLOOD PRES.-LOW(A)
148	AURA DISTORTION	202	BLOOD PRES.-LOW(B)
149	AURA CO-ORDINATION	203	BLOOD SUPPLY-MUSCLES
150	BACH FLOWER AGRIMONY	204	BLOOD(ALTERNATIVE)
151	BACH FLOWER ASPEN	205	BONE CONDITIONS
152	BACH FLOWER BEECH	206	BONE(DETERIORATION)
153	BACH FLOWER CENTAURY	207	BONE (PERIOSTEUM)
		208	BONE (MYELOCYTES)
		209	BONE(ERYTHROBLASTS)
		210	BONE(OSTEOGENETIC)

211	BONE REGENERATION	268	CHOLESTEROL ADRENAL
212	BREASTS(STRUCT.TONE)	269	CHOLINE
213	BRIGHTS DISEASE	270	CIRCULATION
214	BRONCHIAL TUBES	271	COLITIS
215	BRONCHITIS	272	COLLOIDIAL SILVER
216	BRUISE	273	COLON
217	BRUISING REDUCTION	274	COLON WALL (1)
218	BURSITIS	275	COLON WALL (2)
219	CALCIFICATION(A)	276	COLON BLOCKAGE
220	CALCIFICATION(B)	277	COLOR YELLOW LUMIN.
221	CANDIDA ALBICANS	278	COLOR GREEN GOVERN.
222	CAPILLARIES	279	COLOR BLUE REFRIG.
223	CARCINOMAS-USE W/ALL	280	COLOR YELLOW VIT. A
224	CAR. ADENOCARCINOMA	281	COLOR RED/OR VIT. B
225	CAR. ASTRICYTOMA	282	COLOR LEMON VIT. C
226	CAR.CHORION EPITHEL.	283	COLOR VIOLET VIT. D
227	CAR. ENCEPHALOID	284	COLOR SCARLET VIT. E
228	CAR.GLIOMASTOMA MUL	285	COLOR MAGENTA VIT. E
229	CAR.LYMPHATIC LEUKEM	286	COLOR INDIGO VIT. K
230	CAR.HODGKINS LYMPHOM	287	COMMON COLD A
231	CAR.NON HODGKINS LYM	288	COMMON COLD B
232	CAR. MELANOMA	289	CONGESTION INTERSTIT
233	CAR.MONOCYTIC LEUKEM	290	CONGESTIVE HEADACHE
234	CAR.MUCOUS CARCINOMA	291	CONNECTIVE TISSUES
235	CAR. MYELOID SARCOMA	292	CONSTIPATION(A)
236	CAR. NEUROBLASTOMA	293	CONSTIPATION(B)
237	CAR.OSTEOGENIC SARCO	294	CONSTIPATION(C)
238	CAR. RENAL CARCINOMA	295	CONSTRICTION
239	CAR.SPINDLE CELL SAR	296	CONTRACTION-TRISMUS
240	CAR.SCIRRHOUS CARCIN	297	CORTISONE
241	CAR. SPHEROID GLAND	298	COUGH
242	CAR. SYNCYTIOMA	299	CRAMP
243	CAR. CERVICAL	300	CRAMP PAIN(LEGS)CIRC
244	CAR. PROSTATE	301	CYSTITIS(A)BLADDER
245	CARDIO-VASCULAR	302	CYSTITIS(B)
246	CARTILAGE DAM.KNEE	303	DEAFNESS/HEARING
247	CATARRH	304	DEBILITY(A)
248	CATARACTS DISRUPTIVE	305	DEBILITY(B)
249	CATARACTS REPAIR LEN	306	DENTAL GUMS
250	CELL REGENERATION A	307	DENTAL TEETH (A)
251	CELL REGENERATION B	308	DENTAL TEETH (B)
252	CELLS-CLEAN	309	DIAPHRAGM
253	CELLULITIS	310	DISLOCATION
254	CEREBRAL SPINAL FLD.	311	DISPLACEMENT
255	CEREBRAL CORTEX	312	DISSEM.SCLEROSIS(A)
256	CERVIX	313	DISSEM.SCLEROSIS(B)
257	CHAKRAS-BASE	314	DISTENSION
258	CHAKRAS-SACRAL	315	DIVERTICULITIS
259	CHAKRAS-SOLAR PLEXUS	316	DOWNS SYNDROME (A)
260	CHAKRAS-HEART	317	DOWNS SYNDROME (B)
261	CHAKRAS-THROAT	318	DRUGS
262	CHAKRAS-BROW	319	DUPUYTRENS CONTRACT.
263	CHAKRAS-CROWN	320	DUODENUM
264	CHICKEN POX	321	DYSPEPSIA
265	CHOLESTEROL NERVE FB	322	EAR-INNER
266	CHOLESTEROL LIVER	323	EAR-MIDDLE
267	CHOLESTEROL KIDNEY	324	EDEMA

325	ELIMINATE FLUIDS ETC	382	HEAD-BRAIN COR.COL.A
326	ELIMINATE ALL METALS	383	HEAD-BRAIN COR.COL.B
327	ELIMINATE THE CAUSE	384	HEAD-BRAIN COR.COL.C
328	EMPHYSEMA	385	HEAD-CRANIAL OSTEOP.
329	ENDOMETRIOSIS	386	HEAD-CROWN
330	ENDOCRINE-PINEAL	387	HEAD-NOISES
331	ENDO.PINEAL(3RD EYE)	388	HEAD-HYPOTHALAMUS
332	ENDO.-PITUITARY(WHL)	389	HEAD-THALAMUS
333	ENDO.-PITUITARY(EXT)	390	HEARING
334	ENDO.-PITUITARY(POS)	391	HEART
335	ENDOCRINE-THYROID A	392	HEARTBURN
336	ENDOCRINE-THYROID B	393	HEMORRAGE
337	ENDO.THYROID FUNCTN.	394	HEPATITIS
338	ENDO.THYROID HORMONE	395	HERNIA
339	ENDO.-PARATHYROID	396	HERPES ZOSTER(A)
340	ENDOCRINE-THYMUS	397	HERPES ZOSTER(B)
341	ENDO.-ADRENAL(WHOLE)	398	HODGKINS DISEASE
342	ENDO.-ADRENAL(CORTX)	399	HOOKWORMS
343	ENDO.-ADRENAL(MEDUL)	400	HORMONE IMBALANCE
344	ENDOCRINE-PROSTATE	401	HYPERSENSITIVITY
345	ENDO.-COCCYNGEAL	402	HYPER/OVER SENSITIVE
346	ENDOCRINE-GENERAL	403	ICHTHYOSIS
347	ENDOCRINE-PANCREAS	404	IMMUNE SYS.-BLOOD
348	ENERGY	405	IMMUNE SYS.-SPLEEN
349	EPSTEIN BARR VIRUS	406	IMPINGED NERVE
350	ETHERIC BODY	407	INDIGESTION
351	EXCESSIVE WEAKNESS	408	INFLAMMATION (1)
352	EXHAUSTION	409	INFLAMMATION (2)
353	EYE	410	INFLAMMATION (3)
354	EYES AND EYELID	411	INFLUENZA (1)
355	EYES EPITHELIAL CELL	412	INFLUENZA C (2)
356	FALLOPIAN TUBES	413	INFLUENZA (3)
357	FATIGUE	414	INFLUENZA(A)
358	FAT METABOLISM	415	INFLUENZA(B)
359	FATTY TISSUE	416	INJURY
360	FATTY TUMOUR	417	INJURY(MUSCULAR)
361	FATTY TUMOUR(LIPOMA)	418	INJURY(BONE)
362	FIBROIDS (A)	419	INSOMNIA
363	FIBROIDS (B)	420	INTESTINE(SMALL)
364	FIBROMA	421	INTESTINE(LARGE)
365	FIBROMYALGAE	422	JOINTS
366	FIBROSIS	423	KIDNEYS
367	FIBROSITIS	424	KNEE INTEROPEOUS
368	FLATULENCE	425	LASSITUDE
369	FOLLICLE	426	LESION
370	FORGIVENESS	427	LIGAMENT (A)
371	FRACTURE	428	LIGAMENT (B)
372	GALL BLADDER	429	LIPOMA
373	GASTRIC INFLUENCE	430	LIVER
374	GASTRIC ULCER	431	LIVER-CLEAN
375	GASTRITIS	432	LUMBAGO
376	GLANDULAR FEVER	433	LUNGS
377	GLAUCOMA	434	LUPUS VULGARIS (A)
378	GOUT	435	LUPUS VULGARIS (B)
379	HAIR	436	LUPUS ERYTH (A)
380	HEAD-BRAIN	437	LUPUS ERYTH (B)
381	HEAD-BROW	438	LUPUS VULGARIS6/96 C

439	LUPUS ERYTH DISCOID	495	MUSCLE JOINT
440	LUPUS ERYTH SYSTEMIC	496	MUSCLE CELLS
441	LYMPHATIC DUCTS	497	MUSCLE(INVOLUNTARY)
442	LYMPHATIC GLANDS	498	MUSCLE(VOLUNTARY)
443	LYMPHATIC LEUKAEMIA	499	MUSCLE FIBRES
444	LYMPHATICS	500	MUSCULAR(ELIM.PAIN)
445	MALTA FEV/BRUCellos.	501	MUSCLE TONE(OPTIM)
446	ME	502	MUSCLE TONE(STREN.)
447	MEMBRANE MUCOUS	503	MUSCULAR ATROPHY
448	MENIERES DISEASE	504	MYELTIN SHEATH
449	MENOPAUSAL-OESTROGEN	505	MYELOID
450	MENOPAUSAL-PROGEST.	506	NAILS
451	MENOPAUSAL-TESTOST.	507	NASAL CONDITIONS
452	MENOPAUSE (A)	508	NECK
453	MENOPAUSE (B)	509	NERVE COND.(DEPR.) B
454	MENSTRUATION	510	NERVE COND.(DEPR.) A
455	MERIDIAN-LUNGS	511	NERVES
456	MERIDIAN-LG.INTEST.	512	NERV.-ENCEPHALON
457	MERIDIAN-PERICARDIUM	513	NERV.-HYPOTHALAMUS
458	MERIDIAN-3WARM-HIGH	514	NERV.-CARESELLUM
459	MERIDIAN-3WARM-MID.	515	NERV.-SYMPATHETIC
460	MERIDIAN-3WARM-LOW	516	NERV.-PARASYMPATHET.
461	MERIDIAN-HEART	517	NERV.-PARASMP./SYMP.
462	MERIDIAN-SM.INTEST.	518	NERVE TONE
463	MERIDIAN-SPLEEN/PANC	519	NERVE CELL
464	MERIDIAN-LIVER	520	NERVE FIBRE
465	MERIDIAN-STOMACH	521	NEURITIS
466	MERIDIAN-CIRC.SEX	522	NEUROMA
467	MERIDIAN-GALL BLAD.	523	NEURONE AXON NV.FIB.
468	MERIDIAN-KIDNEYS	524	NOSE
469	MERIDIAN-BLADDER	525	OBESITY(CELLU.RED)
470	MERIDIAN-SKIN	526	OPTIC NERVE
471	MERIDIAN-NERV.SYST.	527	OVARIES(A)
472	MERIDIAN-LYMPH.	528	OVARIES(B)
473	MERIDIANS	529	OXYGEN
474	METALS CHROME	530	OXYGEN-TISSUE
475	METALS COPPER	531	OZONE
476	METALS GOLD	532	PARASITES-GENERAL
477	METALS LEAD (POISON)	533	PARKINSONS DISEASE
478	METALS MERCURY	534	PENIS (STRICTURE)
479	METALS MOLYBDENUM	535	PERIOSTEUM
480	METALS NICKEL	536	P.I. PHYSICAL INJURY
481	METALS PALLADIUM	(* Car. Test on later Mk6s)	
482	METALS PLATINUM	537	PHYSICAL INJURY-BONE
483	METALS SILVER	(*Car. Leukaemia on later Mk6s)	
484	METALS TIN	538	PHYSICAL INJURY-MUSC
485	METALS TITANIUM	(Melatonin on later Mk6s)	
486	METALS ZINC	539	PHYSICAL RELAXATIONS
487	MIGRAINE	540	PHLEBITIS
488	MINERALS CALCIUM	541	PLANAR WARTS
489	MINERALS MAGNESIUM	542	PLASMODIUM MALAISE
490	MINERALS POTASSIUM	543	PLURAL
491	MONEY	544	PNEUMONIA CONDITION
492	MONOCYTIC LEUKEMIA*	545	POLIO
(*Motor Neurons on later Mk6		546	POLYPUS
493	MUSCLE	547	PONS VAROLIUS
494	MUSCLE STIMULATION	548	POISON LEAD

549	POISON HEROIN	606	SEROTONIN(5-HYDROXY)
550	POISON HASHISH	607	SHOCK(MENTAL)(ELECT)
551	POISON DETOX.	608	SHOCK (PHYSICAL)
552	POISON TOBACCO	609	SHOULDER JOINT
553	POSITIVE ALT.CURRENT	610	SILICONE
554	PSORIASIS	611	SINUS
555	PSYCH.-ALPHA STATE	612	SKIN DISEASES
556	PSYCH.-CALMNESS	613	SKIN (PSORIASIS)
557	PSYCH.-CHEERFULNESS	614	SKIN(GENERAL/DRY)
558	PSYCH.-CONCENTRATION	615	SMELL
559	PSYCH.-CONFIDENCE	616	SNAKE BITE
560	PSYCH.-CONSTR.DESIRE	617	SPERMATIC CORDS
561	PSYCH.-COURAGE	618	SPINE (DISCS) (A)
562	PSYCH.-ECSTASY	619	SPINE (DISCS) (B)
563	PSYCH.-ENTHUSIASM	620	SPINE LIGAMENTS (A)
564	PSYCH.-FOREBEARANCE	621	SPINE LIGAMENTS (B)
565	PSYCH.-FORCEFULNESS	622	SPINAL NERVES
566	PSYCH.-HAPPINESS	623	SPINE (VERT.JOINT)
567	PSYCH.-HOPEFULNESS	624	SPINE (VERTEBREA)
568	PSYCH.-HYPERTENSION	625	SPINETER
569	PSYCH.INDUSTRY	626	SPLEEN
570	PSYCH.-LOVE	627	SPONDYLITIS
571	PSYCH.-MELANCHOLIA	628	SPRAIN
572	PSYCH.-MEMORY	629	STAPHYLOCOCCUS
573	PSYCH.-MENTAL ALERT.	630	STIFFNESS (MUSCLE)
574	PSYCH.-MENTAL DEPRS.	631	STIMULATION-GENERAL
575	PSYCH.-MENTAL DRIVE	632	STOMACH
576	PSYCH.-MENTAL ENERGY	633	STOP SMOKING
577	PSYCH.-MENTAL RELAXN	634	STRAIN(MUSCULAR)
578	PSYCH.-MENTAL VITALY	635	STREPTOCOCCUS
579	PSYCH.-NERV.DEPRESS.	636	STRESS RELIEF
580	PSYCH.-PERSONAL MAG.	637	STROKE
581	PSYCH.-SHOCK(MENTAL)	638	SUBLUXATION
582	PSYCH.-SEX DES. MALE	639	SWEAT GLANDS
583	PSYCH.-SEX DES. FEM.	640	SWELLING(REDUCTION)
584	PSYCH.-SEX FUNT.MALE	641	SYNOVIAL MEMBRANE
585	PSYCH.-SEX FUNT.FEM.	642	SYNOVIAL BURSAE
586	PSYCH.-SKIN DISEASES	643	TASTE
587	PSYCH.-STAMINA	644	TEETH
588	PSYCH.-STRENGTH	645	TENDON
589	PSYCH.-SUICIDAL	646	TESTICLES(OCHAITIS)
590	PSYCH.-GEN.TENSION	647	THROAT(SORE THROAT)
591	PSYCH.-TENSION	648	TINITUS
592	PSYCH.-VITAL LIFE FO	649	TISSUE DAMAGE(TRAU.)
593	RHEUMATISM	650	TMJ
594	RHEUMATISM-ACUTE	651	TONSILS
595	RHEUMATISM-CHRONIC	652	TORTICOLLIS SPASM.
596	RHEUMATISM-GONORRHEA	653	TOXIC METALS
597	RHEUMATISM-MUSCULAR	654	TOXICITY (A)
598	RHEUMATISM-HANDS	655	TOXICITY (B)
599	SCAR TISSUE-(A)	656	TOXICITY (C)
600	SCAR TISSUE-(B)	657	TOXICITY (D)
601	SCAR TISSUE-REGEN.PO	658	TOXICITY (E)
602	SCIATICA (A)	659	TOXICITY (F)
603	SCIATICA (B)	660	TRIGERMINAL NERVE
604	SCIATICA (ALT.COND.)	661	ULCERATION (MOUTH)
605	SEPSIS	662	UTERUS (PROLAPSE)

663 VAGINA
664 VARICOSE ULCER
665 VARICOSE VEINS
666 VEINS (PRE-VARICOSE)
667 VERTEBRAE(DAM./TRAU)
(*Menopause Oestrogen B on later Mk6s)
668 VERTEBRAL JOINTS
(*Menopause Progest. B on later Mk6s)
669 VERTIGO
670 VITALITY(A)
671 VITALITY(B)
672 VITALITY (STABILISE)
673 VITAMIN A1
674 VITAMIN A2
675 VITAMIN B1
676 VITAMIN B2
677 VITAMIN B3
678 VITAMIN B4
679 VITAMIN B5
680 VITAMIN B6
681 VITAMIN B12
682 VITAMIN C
683 VITAMIN D1
684 VITAMIN D2
685 VITAMIN D3
686 VITAMIN E
687 VITAMIN F
688 VITAMIN G
689 VITAMIN H
690 VITAMIN K1
691 VITAMIN K2
692 VOCAL CORDS(STR.DAM)
693 VOICE(LARYNGITIS)
694 WATER-ANTI BACTERIA
695 CALIBRATION-10Hz
(*Psych. Delta Waves on later Mk6s)
696 CALIBRATION-100Hz
(*Psych. Theta Waves on later Mk6s)
697 CALIBRATION-1000Hz
(*Radiation Exposure on later Mk6s)
698 CALIBRATION-9999Hz
(Calibration 10Hz on later Mk6s)